



The BD Integra™ Syringe is available in an array of sizes for all SubQ and IM injections.

Available Fall/Winter 2002

BD Integra 3ml Syringe w/ Detachable Needle

Catalog #	Description	Gauge/Needle	Packaging
305269	3ml	25G, 5/8"	100/Box – 400/Case
305270	3ml	25G, 1"	100/Box – 400/Case
305271	3ml	23G, 1"	100/Box – 400/Case
305272	3ml	22G, 1-1/2"	100/Box – 400/Case
305273	3ml	21G, 1"	100/Box – 400/Case
305274	3ml	21G, 1-1/2"	100/Box – 400/Case

BD Integra 1ml Syringe w/ Perm. Attached Needle

Catalog #	Description	Gauge/Needle	Packaging
305282	Insulin	29G, 1/2"	100/Box – 800/Case
305281	Insulin	28G, 1/2"	100/Box – 800/Case
305298	TB	27G, 1/2"	100/Box – 800/Case
305279	TB	26G, 3/8"	100/Box – 800/Case
305278	TB	25G, 1"	100/Box – 800/Case
305277	TB	25G, 5/8"	100/Box – 800/Case

BD Integra Needles, 3ml Syringe

Catalog #	Description	Gauge/Needle	Packaging
305310	Retracting Ndl.	25G, 5/8"	100/Box – 1000/Case
305311	Retracting Ndl.	25G, 1"	100/Box – 1000/Case
305312	Retracting Ndl.	23G, 1"	100/Box – 1000/Case
305313	Retracting Ndl.	22G, 1-1/2"	100/Box – 1000/Case
305275 ¹	Blunt Fill Needle	18G, 1-1/2"	100/Box – 1000/Case
305283	3ml Syringe Only		100/Box – 400/Case

¹ 18G Blunt Fill Needle has a non-retracting needle

For more information call
1.888.237.2762



BD, BD Logo, and all other trademarks are property of Becton, Dickinson and Company.
© 2002 BD.

BD Medical Systems
1 Becton Drive
Franklin Lakes, NJ 07417
www.bd.com



Indispensable to
human health

Shaping the Future of Safety

BD Integra™ Syringe

with Retracting BD PrecisionGlide™ Needle

Step-by-Step Safety

1

Aspiration

- Draw up medication in accordance with established protocol.

2

Injection

- Administer medication in accordance with established protocols, making sure that all medication has been dispensed.

3

Activation

- Push plunger rod to activate the safety mechanism. Depress until two “clicks” are sensed.

* Activation can occur before or after the needle is withdrawn from the patient.

4

Needle Retraction

- The needle will disappear into the syringe and the plunger will be recessed into the barrel when activation is complete.

Discard immediately in nearest sharps container.

The BD Integra™ Syringe Advantage:

Detachable Needle¹



- Allows Clinicians to change needles for aspiration and administration
- Uses Tru-Lok Technology which provides:
 - Low waste space
 - Secure connection

Activates before or after needle is withdrawn from patient

- Clinician makes the decision
- Complete dosing ensured prior to activation


Uses BD PrecisionGlide™ “Sharpest” Needle Technology²

- Maintains its sharpness
- Improves patient comfort


¹Detachable Needle only Available on 3ml syringe with Needle.
²Report on the Evaluation of the Improved BD PrecisionGlide™ Needle, “The World’s Sharpest Needle”.

The BD Integra™ Syringe combines innovative technology with the BD quality and comfort you trust.

PRIOR TO ACTIVATION

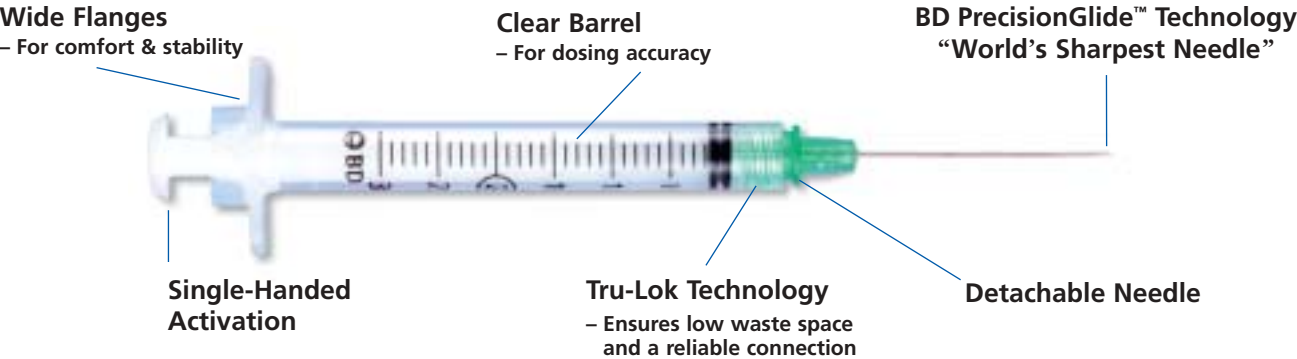


AFTER ACTIVATION

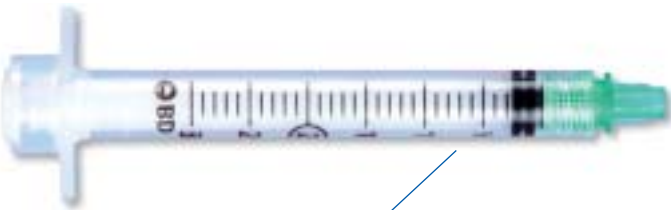

Needle Retracts into the Syringe!
Eliminates exposure to the needle after activation

Because of its Tru-Lok design, the BD Integra™ Syringe has one of the lowest waste space volumes.

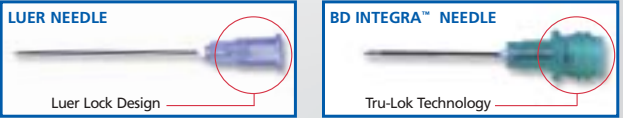

Tru-Lok is a unique, non-luer thread design specific to the BD Integra™ Syringe and its Detachable Needles. With Tru-Lok Technology, the BD Integra™ Syringe and Needle achieves .026ml of waste space.

Low waste space means reduced medication waste and medication costs for a healthcare facility. The example below illustrates the variation in cost and number of doses obtained from a multi-dose vial when using a BD Integra™ Syringe, BD Conventional Syringe or Brand A Retracting Syringe.

HEPATITIS B CHART			
SYRINGE TYPE	DRUG	# OF DOSES PER VIAL	COST PER DOSE
BD Integra™ Syringe	Hepatitis B	3.18	\$65.94
BD Conventional Syringe	Hepatitis B	3.00	\$69.82
Brand A Retracting Syringe	Hepatitis B	2.70	\$77.57

The BD Integra™ Syringe provides a medication cost savings!

Example:
 $\$69.82 - \$65.94 = \$3.88$ (a dose) x 1000 (doses) = **SAVE! \$3,880.**