

DON'T FORGET...

Aseptic technique and proper skin preparation of the site are essential.

Discard BD Safety-Lok™ Syringe in the nearest sharps collector. Follow your institution's policy for safe disposal of all medical waste.

Safety-Lok™ Offerings		
CATALOG #	NEEDLE	PACKAGING
1ml Insulin Syringe with Permanently Attached Needle		
329464	29G 1/2" U-100 Ultra-Fine	100/Box - 500/Case
1ml Tuberculin Syringe with Permanently Attached Needle		
305553	27G 1/2"	100/Box - 500/Case
305554	25G 5/8"	100/Box - 500/Case
3ml Syringe with Detachable Needle		
309592	25G 5/8"	100/Box - 800/Case
309594	23G 1"	100/Box - 800/Case
309593	22G 1-1/2"	100/Box - 800/Case
309596	22G 1"	100/Box - 800/Case
309595	21G 1-1/2"	100/Box - 800/Case
5ml Syringe with Detachable Needle		
305561	21G 1-1/2"	50/Box - 400/Case
10ml Syringe with Detachable Needle		
305564	21G 1-1/2"	50/Box - 400/Case
Syringes Only, Luer-Lok® Tip		
309606	3 ml Syringe Only	100/Box - 800/Case
305558	5 ml Syringe Only	50/Box - 400/Case
305559	10 ml Syringe Only	50/Box - 400/Case

1.800.817.1702



BD Medical Systems
1 Becton Drive
Franklin Lakes, NJ 07417-1884
www.bd.com

Becton Dickinson Canada Inc.
2771 Bristol Circle
Oakville, Ontario L6H 6R5
www.bd.com

SafetyGlide and Safety-Lok are trademarks of Becton, Dickinson and Company. Luer-Lok is a registered trademark of Becton, Dickinson and Company.
BD 11/99-0403
©1999 Becton, Dickinson and Company



The World's First Safety Syringe

BD Safety-Lok™
Sliding Sleeve Syringe

Step-by-Step Safety

1 Assembly
To assemble needle and syringe, hold syringe by flanges (wings) and twist needle until firmly seated.

2 Preparation
During aspiration, press index finger against flanges to prevent sleeve movement.

3 Post-Injection
After injection, grasp sleeve firmly and twist flanges to loosen the sleeve.

4 Lock Fully retract needle into sleeve until it locks in the protected position. When the Safe Lock indicator green band fully covers the red band and an audible click is heard, the sleeve is locked into position.

Establishing the first generation.

BD's record of innovation spans a century. Time after time, BD has led the way with new products and technologies, especially in healthcare worker safety. More than 10 years ago, BD pioneered safety injection with the introduction of the Safety-Lok™ Syringe.



Today, BD Safety-Lok™ is the best-selling safety hypodermic device on the market and has allowed many healthcare institutions and practices to provide safety-engineered protection to their clinicians.



Prior to Activation


World's Sharpest Needle
for greater patient comfort

Bold Scale Markings
facilitate accurate dose measurement

Crystal Clear Barrel and Sleeve
make syringe contents easier to see

Large Flanges and Tapered Plunger Rod
increase control for easier aspiration and injection

Efficient Blister Packaging
minimizes storage space and waste disposal costs

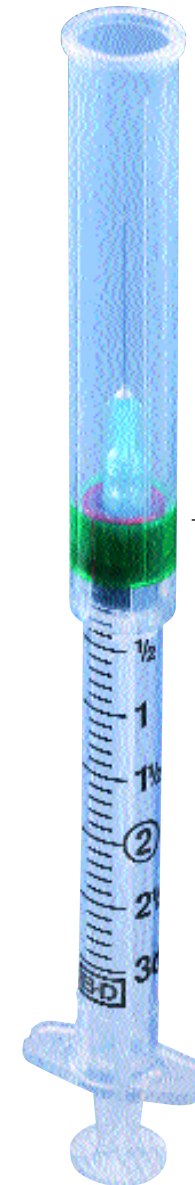

After Activation

Key to Safety


BD Safety-Lok™ is the only safety-engineered syringe of its type to offer clear visual confirmation of the lock position.

Visual Lock Indicator

Lock confirmed when green band fully covers red collar.



Complete Lock

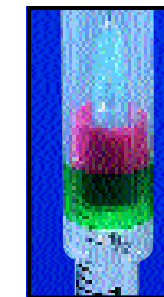

Incomplete Lock

Pioneering the next generation.

At BD, the next generation of safety injection devices has arrived – the BD SafetyGlide™ Shielding Hypodermic Needle. Developed based on extensive clinical feedback from hundreds of interviews, focus groups and clinical immersions, BD SafetyGlide™ is setting a new standard in safety injection.*

BD SafetyGlide™ offers single-handed, immediate activation with minimal change in technique. Since the BD SafetyGlide™ detachable needle fits onto conventional syringes, including Luer-Lok® pre-fills, stocking requirements are minimized.

* Current high demand limits availability until Summer 2000.



BD. First in Safety.